

REDUCE THE RISK OF FALLS & INJURIES

LOW PROFILE, NON-SLIP MATERIAL

PORTABLE

EASY TO USE

EASY TO MOVE FROM ROOM TO ROOM

WATERPROOF

WASHABLE

INTEGRATES TO MOST NURSE CALL SYSTEMS


PLUG & PLAY OPERATION

DURABLE & RELIABLE

COST EFFECTIVE

SMALL PRESSURE MAT (380mm X 255mm)

Chair Sensor Pad

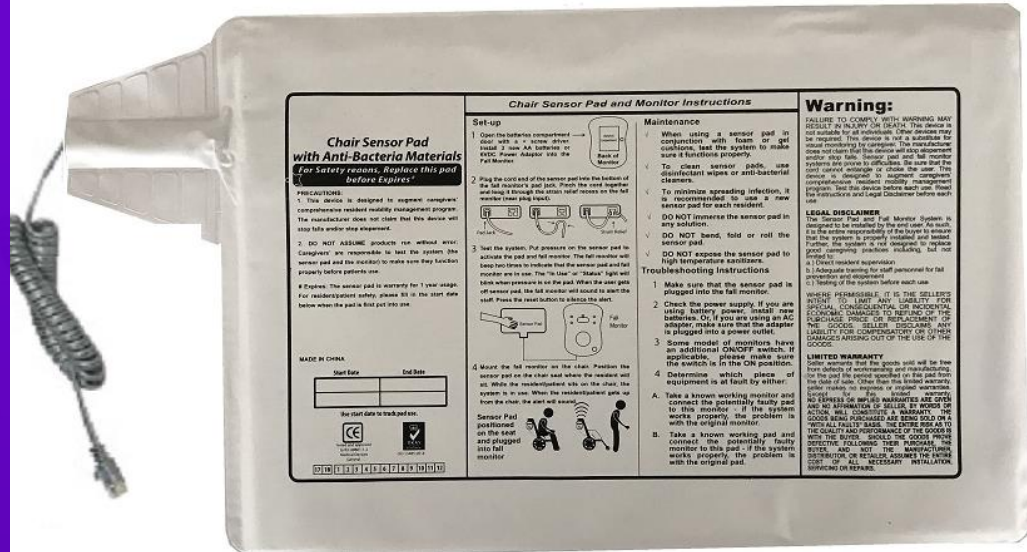

MED-NC-C365-NC10 (with RJ11)
MED-NC-C365-JS63 (with Jack 6.35mm)
(TGA Registration: 223902)

Indigo Care provides a pressure sensitive Premium Chair Sensor Pad that is designed to immediately alert staff should a resident leave their chair. This is most useful to residents who are at risk of falling out of their chair, or those who require assistance in arising from the chair. The pad is placed in the centre of the chair, directly underneath the resident's bottom. If a chair cushion is used, then the Sensor Pad must be placed on top of the cushion. For a normally open call system, a monitor is connected to the chair sensor pad and then, with a Y cable (via a monitor). Once the resident removes themselves from the chair pad and tries to leave their chair, the sensor mat will activate the monitor's alarm and alert the staff. Set-up instructions are printed on the actual chair pad.

KEY SPECIFICATIONS:

- **Dimensions: 380mm X 255mm**
- **Connection: RJ11 on 500mm curly cord (normally open) or Jack 6.35mm plug (normally close)**
- **Silver Printing, Resistance <1k Ohm**
- **White PVC, Thickness 0.20mm**
- **IP Rating : 65**
- **Warranty: 12 months**

PRODUCT ACCESSORIES:

- **RJ11 requires a Y Cable**
- **Jack 6.35 requires straight 6.3mm or 3.5mm cable**
- **Anti-Tamper Monitor (use 4 x AA batteries-not included)**
- **6VDC Adaptor (optional)**

Horizon Electronic Components T/as INDIGO CARE
 2/71 Rushdale Street, Knoxfield, VIC, 3180
 PH +61 3 9763-2988 FAX +61 3 9763-3811
 www.indigocare.com.au
 ABN 36007302801



ORDER INFORMATION

Fax or email your purchase orders to
03-9763-3811
 Or email sales@indigocare.com.au

